

INSTRUCTIONS
FOR
MODEL NO. HT-300/500 Ver: 3.0
HARDNESS TESTER

Manual Release: 3.0

KEY INTERNATIONAL, INC.
480 ROUTE 9
ENGLISHTOWN, NJ 07726
PHONE: (732) 536-9700 FAX: (732) 972-2630

Publication Record

Release No.	Date	Description
1.0	01/30/1999	Original Printing
2.0	05/19/1999	Added the instructions for the Vari-Flex jaws to the manual. Fixed an error with regard to the clear key. Changed the description of the Start key operation.
3.0	09/27/1999	Removed the old keyboard part number from "Spare Parts List" and replaced with the new keyboard part number. Added the Main Board part number to list.

IMPORTANT NOTICE

THIS MODEL HT-300/500 TABLET HARDNESS TESTER HAS BEEN EQUIPPED WITH NEWLY DESIGNED ELECTRONICS. THESE CHANGES ARE PRESENT IN ALL MODELS AFTER SERIAL NUMBER I-4500/B. THE FUNCTION OF THE MACHINE IS UNCHANGED. HOWEVER, DUE TO THE NEW ELECTRONICS THE METHOD OF CALIBRATION AND UNITS SELECTION HAS BEEN CHANGED. THE CHANGES ARE DESCRIBED WITHIN THIS MANUAL.

KEY INTERNATIONAL, INC.
480 ROUTE 9
ENGLISHTOWN, NJ 07726
PHONE: (732) 536-9700 FAX: (732) 972-2630

TABLE OF CONTENTS

1. General Description	5
1.1 Introduction.....	5
1.2 Specifications.....	5
1.3 Method Of Operation.....	6
2. Putting Into Operation.....	7
2.1 Wiring	7
2.2 Inserting The Chart Paper	8
3. Operating keyboard.....	9
3.1 Operating Controls.....	9
4. Test Program.....	13
4.1 Units Selection	13
4.2 Zero Offset	13
5. Calibration Program.....	15
5.1 Validation.....	15
5.2 Calibration.....	16
6. Cleaning Procedure.....	17
7. Vari-Flex Jaw operating instructions	18
7.1 General Description	19
7.2 Installation.....	19
7.3 Setup	20
7.4 Operation.....	20
8. Technical Information.....	22
8.1 Units Conversion List	22
8.2 RS-232 Serial Port	22
8.2.1 Interface	22
8.2.2 Pin Assignments.....	22
8.2.3 Data Output.....	22
8.2.4 Data Input.....	23
9. Optional and Spare parts list.....	24

1. GENERAL DESCRIPTION


1.1 Introduction

The hardness tester model HT-300/500 is a high performance instrument, which can accurately measure the tensile strength of tablets, ampoules, food and other articles, in either single or series operation. The tensile strength can be printed out in Kilo Ponds (KP), Strong Cobbs (SC) or Newtons (N), according to the operator's wish.

1.2 Specifications

SAMPLE		4 - 40 mm
STROKE		4 - 40 mm (externally adjustable)
CHART SPEED		3 mm/sec.
FORCE INCREASE		20 Newton's (N) per second
MEASURING RANGE	HT-300 HT-500	5 - 300 Newtons (N) 8 - 500 Newtons (N)
MEASURED VALUE		Displayed on a 3-digit LED display and printed on the built-in printer
RESOLUTION		1 Newton (N)
ACCURACY		+/- 1% Full Scale (FS)
MEASURING JAWS		Interchangeable for different sample forms (i.e. tablets, ampoules etc.)
DIMENSIONS		530 mm x 240 mm x 180 mm (W x D x H)
MAINS CONNECTION		100 - 120 Volt, 50/60 Hz. switch selectable to 200 - 240 Volt, 50/60 Hz.
I/O INTERFACE		RS-232 Serial

1.3 Method Of Operation

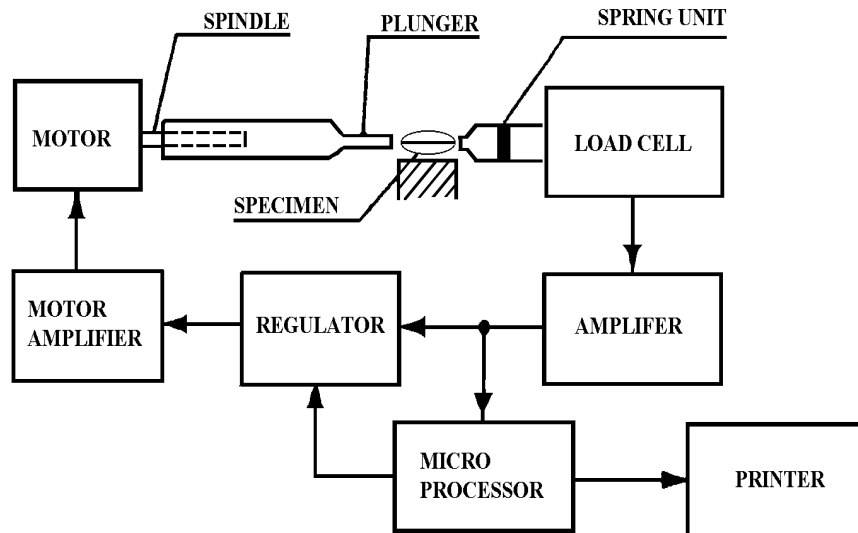


After pressing the START key, the zero-offset value of the load cell is determined and stored in the microprocessor. If the value is within the permitted range, the motor, controlled by the regulator unit and the motor amplifier, causes the plunger to move at a constant speed towards the specimen. As soon as the plunger touches the specimen and produces a force on the load cell, the microprocessor switches the regulator on so that a linearly increasing force is produced via the control circuit (load cell - regulator - motor - specimen.) The spring unit in front of the load cell provides a smooth transition from the constant velocity phase of the plunger to the linear power increase phase. A decreasing signal from the load cell indicates to the microprocessor that the break has occurred. The microprocessor then subtracts the zero-offset value stored at the beginning of the test from this value and prints out the calculated value in the desired units (Newtons, Kilo Ponds or Strong Cobbs).

2. PUTTING INTO OPERATION

Before putting the instrument into operation, check the voltage setting (220-Volts/115 Volts) and make sure that the operating voltage corresponds to the local mains supply.

2.1 Wiring

- Connect the signal cables to any additional units (i.e. computer)
- Connect the power supply cable to the local mains supply

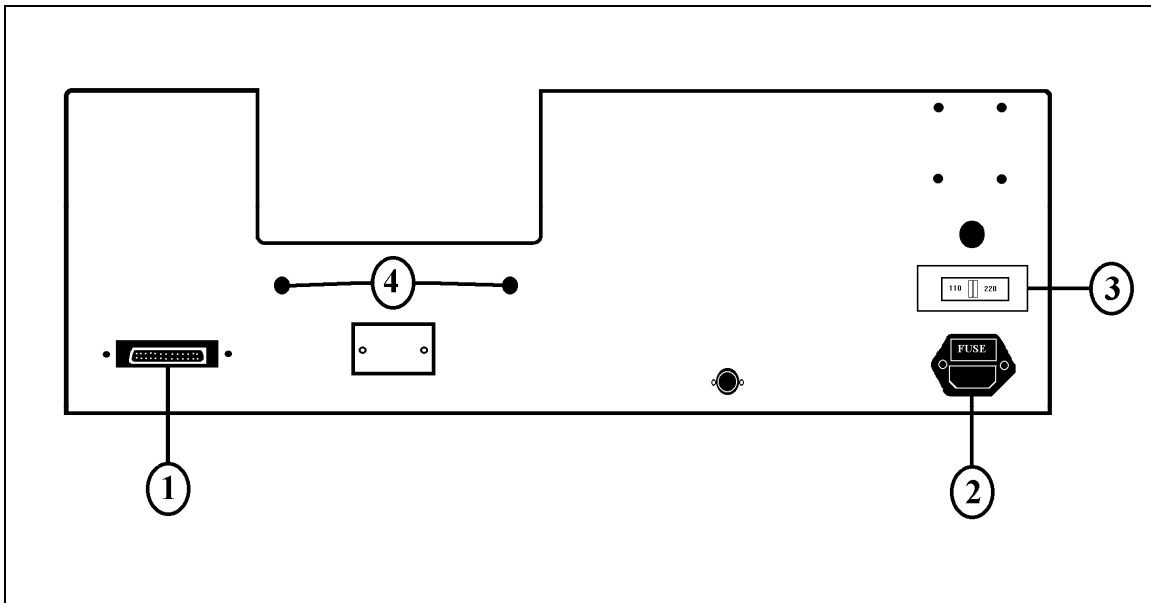

1. I/O Interface 25 pin, Sub-D male connector (RS-232 serial interface)
2. Mains input with fuse (250 mA for 220 Volt operation, 500 mA for 115 Volt operation)
3. Voltage selector switch - 220 volts/115 volts (make sure the correct fuse is installed in the mains connector)
4. Waste container mounting studs

2.2 Inserting The Chart Paper


NOTE: THE INSTRUMENT SHOULD NEVER BE OPERATED WITHOUT PAPER.

To insert the paper, first open the paper door, then remove the paper spindle. Next, slide a roll of paper, with the dim side up, onto the paper spindle and replace the paper spindle in the machine. Now feed the paper into the slot on the black plate of the printer with the dim side facing out. While turning of the paper transport wheel (11) guide the paper through the printer.

When advancing the paper through the printer, do not pull on the paper, use the transport wheel. Pulling on the paper could damage the print head of the printer or cause the paper to jam. If the paper should jam and need to be removed from the printer, be sure to press the paper release button, on the left side of the printer (10), before pulling the paper backward through the printer.

3. OPERATING KEYBOARD

3.1 Operating Controls

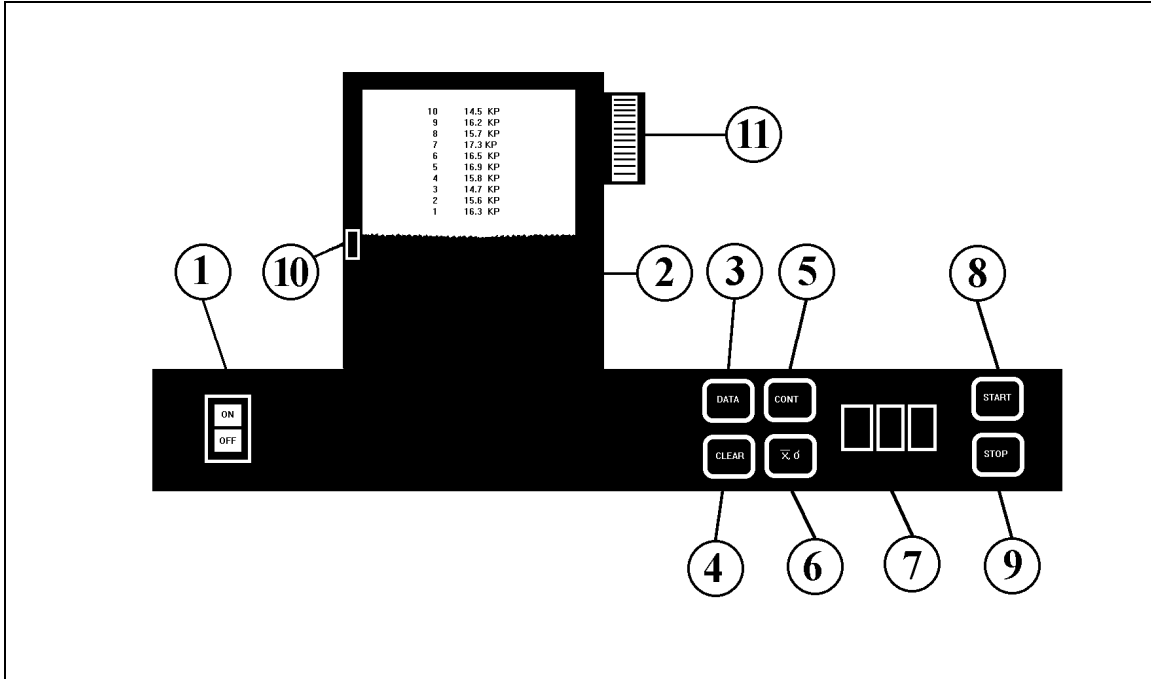

(1) **POWER**

The POWER button (1) is used to switch the instrument's main power on and off.

(2) **PRINTER**

If the printer should jam or the paper tear, wait until the printout is terminated and the printing paper has stopped moving. Move the paper release lever (10) in the direction of the arrow. This will release the paper tension bar and allow the paper to be removed from the bottom of the printer (2). The release lever (10) must not be touched during printer operation. Failure to use the release lever when removing paper from the printer may cause serious damage to the printer.

When advancing the paper through the printer do not pull on the paper. Pulling on the paper may damage the print head or cause a paper jam. Use the thumb wheel to advance the paper through the printer.

Printing paper: Part Number: 15-156-0053
Per box (4 rolls per box)
Per case/36 boxes (144 rolls)

Call for current pricing: Key International, Inc.
480 Route 9
Englishtown, NJ 07726
Phone: (732) 536-9700 Fax: (732) 972-2630

(3) DATA

When the DATA key (3) is activated (LED on), the individual hardness reading of the specimen will be printed on the printout. If switched off the individual hardness readings will not be printed. The DATA key is disabled during the normal testing cycle and will not function until the end of the testing cycle.

(4) CLEAR

If the X, σ key (6) is turned on (LED on), the last measured value can be deleted from the statistical calculation by pushing the CLEAR key (4) (i.e. faulty measurement). The display will show the value that has been cleared and the printer repeats the last line and prints an " = " sign behind the number and result to indicate that the result was cleared.

The CLEAR key (4) is enabled after the break of the specimen and will remain active until the start of the next testing sequence. The specimen that has been cleared will not be included in the statistical printout.

(5) CONT

The CONT key (5) is switched on (LED on) for continuous operation.

Display: "-" no continuous operation.
"1 - 9" Automatic start of next sample occurs from 1 to 9 seconds after the break of the previous specimen.

To change the delay, hold the CONT key (5). By keeping the key pressed the DISPLAY (7) will cycle through the continuous settings. When the setting required is displayed, release the CONT key (5) and the new setting will be in effect.

(6) X, σ

If statistics printout is desired, press the X, σ key (6) before the first measurement. The printer will print "00000000000000000000" and the

LED will turn on. Test the specimens as normal. After the last specimen is tested press the X, σ key (6) again and the statistical report will be printed.

If the statistics report is selected, by pressing the X, σ key (6) (LED on) and the DATA key (3) is turned off (LED off) the statistical report will still print, however, the individual specimen readings will not be printed. If both the X, σ key (6) and the DATA (3) are turned on (both LEDs on), then both the individual specimen values and the statistical report will be printed.

The statistical report includes:

- NR - Number of samples tested
- XM - Mean Value = Sum of Result Values / Number of Samples = $\Sigma X / NR$
- SD - Standard Deviation = $\sqrt{(\Sigma X^2 - NR * XM^2) / (NR - 1)}$
- SREL - Relative Standard Deviation = $SD / XM * 100\%$
- XMIN - Lowest individual hardness value
- XMAX - Highest individual hardness value

The instrument is ready to measure the next specimen as soon as the printout is completed. If the printer has a paper jam or cannot print, the instrument will lock up and no testing can be performed until the problem is solved.

To switch the statistics report mode off, press and hold the CLEAR key (4). While holding the CLEAR key (4), press the X, σ key (6), then release both keys and the LED will turn off.

(7) DISPLAY

The display shows the result of the last measured sample.

In addition the following states are displayed:

"0" after start up.

value the value cleared after a "CLEAR" command.

"-" 1 - 9 on selection of "CONT" operation.

"_" during the measuring procedure.

(8) START

CONT off: By pressing the START key (8), the single measuring procedure is started.

CONT on: By pressing the START key (8), a series of measuring procedures is started. Testing for each progressive specimen after the first will start after the number of seconds set using the CONT key (5) has elapsed.

(9) STOP

By pressing of the STOP key (9), the testing procedure can be interrupted at any time.

(10) PAPER RELEASE

If you need to remove the paper from the printer for any reason, it is very important to push the paper release in the direction of the arrow before pulling the paper backwards through the printer. Failure to do this before removing the paper from the printer may cause damage to the internal mechanism of the printer.

(11) PAPER ADVANCE WHEEL

This wheel is used to advance the paper forward through the printer. Do not pull on the paper to advance it through the printer, as this may damage the print head of the printer or cause a paper jam.

4. TEST PROGRAM

NOTE: BEFORE STARTING THE TEST PROGRAM MAKE SURE THAT THE JAW-STROKE ADJUSTING KNOB HAS BEEN TURNED ALL THE WAY TOWARD THE LARGE CIRCLE. ONCE THIS IS DONE PRESS THE START KEY AND THEN PRESS THE STOP KEY TO CAUSE THE MOVING JAW TO RETURN TO THE NEUTRAL POSITION.

In order to select or change the units in which the tester measures, it is necessary to access the TEST PROGRAM. To enter into the TEST PROGRAM, the movable jaw must be in the home (neutral) position. First, press and hold the STOP key (9) and without releasing it, press the START key (8). The message "PTB-301 XX.XX/XX" will be printed on the printer. Next, release both keys. This number is the version number of the software programmed into the testers hardware and is needed for service and when asking questions about the unit. When this messages is printed, the unit is in TEST PROGRAM mode and all the keys on the keyboard assume new functions, as listed below:

STOP	-	Exits the TEST PROGRAM.
START	-	Units selection.
CONT	-	Display segment test procedure.
X, σ	-	Enters CALIBRATION mode.
DATA	-	Not Used.
CLEAR	-	Prints ZERO OFFSET value.

4.1 Units Selection

After entering the TEST PROGRAM, as described above, press and hold the START key. This will cause the display to cycle through a series of numbers, from zero (0) to three (3). Releasing the START key causes the number display to be selected. After selection, the unit's selected will be printed out on the printer. Each of these numbers represents one of the units of measurement, as listed below:

0	-	Newtons (N).
1	-	Kiloponds (KP).
2	-	Strong Cobbs (SC).
3	-	Digits (Used only for testing purposes).

4.2 Zero Offset

This is the internal value at which the load cell is set to indicate zero (neutral) pressure. This value must be between ten (10) and twenty (20). To verify that the ZERO OFFSET value is correct, press the CLEAR key (4) while the unit is in the TEST PROGRAM mode and the value will be printed on the printer. If the value is not within the specified range, the unit will not function properly and it is necessary to contact the Service

Key International, Inc. • Englishtown, NJ 07726

Department at Key International, Inc.

5. CALIBRATION PROGRAM

The unit needs to be calibrated at least two times a year. Before proceeding with the calibration or validation of the unit it is necessary to have a class F certified weight kit. This kit contains a range of weights, from 5 N to a combined total of 98.066 N that are required for calibration and validation of the unit. The weight kit needs to be re-certified each year and is available from Key International, Inc.

As the certification of the calibration kit can be costly, the unit can be returned to Key International, Inc. for service and calibration. Returning the unit to Key International, Inc. for calibration by a qualified service technician assures that all the functions of the unit are operating correctly.

In order to calibrate the unit, it is necessary to first enter the TEST PROGRAM and then press the X, σ key. After the X, σ key has been pressed, the unit will print a list of the new keyboard functions and activate the CALIBRATION PROGRAM. The new functions are listed below:

STOP	-	Exits the CALIBRATION PROGRAM.
START	-	Validation of load cell.
X, σ	-	Print stored calibration settings in digits.
CONT	-	Prints load cell reading in digits.
CLEAR	-	Activate calibration zero reading. ** SEE NOTE
DATA	-	Activate calibration weight reading. ** SEE NOTE

NOTE: THESE KEYS ALTER THE INTERNAL CALIBRATION READINGS OF THE UNIT AND SHOULD NOT BE USED UNLESS THE VALIDATION READINGS OF THE UNIT ARE INCORRECT. THE PROCESS OF CALIBRATION REQUIRES THE USE OF CERTIFIED WEIGHTS AS DESCRIBED ABOVE.

Once the CALIBRATION PROGRAM has been started, remove the moving jaw by pulling it straight out of the socket. Next, remove the pedestal by unscrewing it from the base of the machine. Now remove the small jaw by pulling it out of the socket. Finally, turn the unit on its side so the black jaw stroke adjusting knob is pointing upward.

5.1 Validation

While in the CALIBRATION PROGRAM as described above it is possible to validate or test the accuracy of the system calibration. It is recommended that a range of weights be tested. The weights should include 0 N, 5 N, 10 N, 20 N, 30 N, 40 N, 50 N, 60 N, 70 N, 80 N and 98 N. Each weight or group of weights is placed on the load cell and the operator then presses the START key. This will cause the printer to print "VALIDATE = XXX N", where "XXX" is the value of the weight tested, in Newtons (N). If any of the values tested are not correct, the unit needs to be re-calibrated.

5.2 Calibration

While in the CALIBRATION PROGRAM, it is necessary to first press the X, σ key. This will print the calibration data stored in the internal circuitry of the unit on the built-in printer. The internal circuitry of the tester measures all values in a unit known as digits. A digit is equal to approximately .348 N. The first line of the printout will show the number of digits for the 0 KP reading of the load cell and the second line will show the number of digits for the calibration weight value of 10 KP (98.066 N). This function should be performed before continuing with the calibration procedure.

Next, with no weight on the load cell, press the CLEAR key. This will cause the internal circuitry to store the new 0 KP value.

Now place the combined weight of 98.066 N on the load cell and press the DATA key. This will cause the internal circuitry to store the new 10 KP (98.066 N) value. The unit has now been calibrated. It is now recommended that the X, σ key be pressed and the new values be compared to the original ones printed at the beginning of the calibration procedure.


When the calibration procedure is finished, the unit must be validated again by following the instructions in the Validation section of this manual.

6. CLEANING PROCEDURE

NOTE: NEVER USE AGGRESSIVE CHEMICAL SOLVENTS THAT MIGHT FOG THE PLASTIC PARTS WHEN CLEANING THE UNIT.

- (1) Before starting the cleaning procedure, make sure the unit is turned off. Next, turn the jaw stroke adjusting knob, on the left side of the machine, completely toward the large circle. Now turn the unit on. If the jaw does not move to the full open position press the START key. As soon as the jaw begins to move forward press the STOP key and the jaws should move to the full open position. Once the jaws are in the full open position, turn the unit off.
- (2) Disconnect the power supply cord. **TO AVOID POSSIBLE ELECTRIC SHOCK DO NOT ATTEMPT TO CLEAN THE UNIT WHILE IT IS STILL CONNECTED TO THE POWER SUPPLY.**
- (3) Remove the long jaw by pulling it out of its socket. Remove the pedestal by unscrewing it from the base of the machine. Remove the short jaw by pulling it out of its socket. If tablet debris has been sticking to the jaws during testing, it may be necessary to polish the contact faces of the jaws.
- (4) When cleaning the unit, it is recommended that isopropyl alcohol be used. **DO NOT ALLOW THE CLEANER TO ENTER THE PRINTER OR REACH THE INTERNAL ELECTRONICS AS THIS MAY CAUSE DAMAGE TO THE UNIT.**
- (5) After the unit has been cleaned, replace the small jaw in its socket. Then screw the pedestal back into the base of the machine and replace the long jaw in its socket.

7. VARI-FLEX JAW OPERATING INSTRUCTIONS

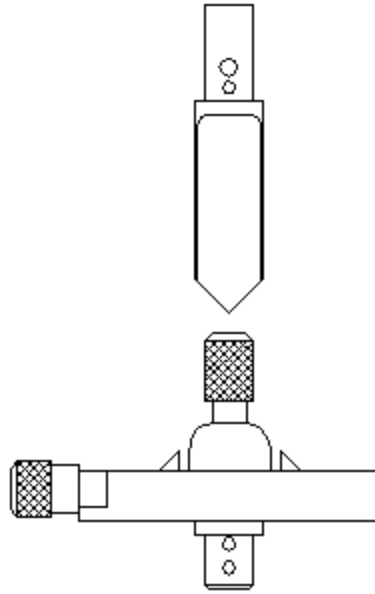


7.1 General Description

The Vari-Flex jaws have been specifically designed for the Key International, Inc. Model HT-300/500 hardness testers. They are a set of two special purpose jaws that replace the standard jaws and pedestal arrangement on the unit. The Vari-Flex jaws provide the user with the ability to break tablets or other products using the "flexure" method of hardness testing. They provide the operator with an adjustable, lead-screw type assembly that allows the jaws to be used on a broad range of product sizes and shapes.

7.2 Installation

To install the Vari-Flex jaws for use in the Model HT-300/500 hardness tester, it is first necessary to remove the standard jaws and pedestal from the unit. This is done by pulling the long jaw straight out of its socket, then unscrewing the pedestal from the base of the machine using a counter-clockwise rotation and, finally, pulling the short jaw straight out of its socket. Once these are removed, it is possible to slide the adjustable tablet holding jaw, hereafter referred to as "Jaw 2", into the vacated socket of the short or stationary jaw. Next, slide the pointed jaw, hereafter referred to as "Jaw 1", into the socket vacated by the long or moving jaw. Once both Vari-Flex jaws have been installed it is possible to proceed to the setup section.


Jaw 1 - Moving jaw (Motor side)

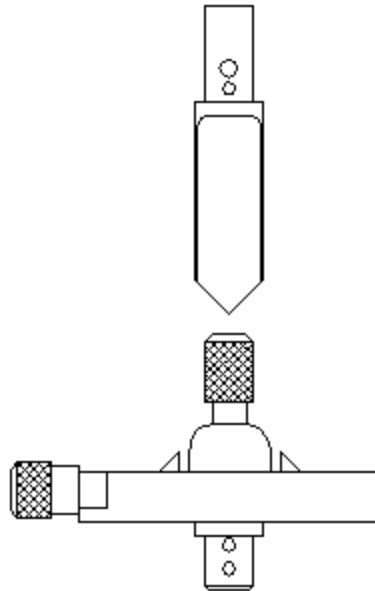


Jaw 2 - Stationary jaw (Load Cell side)

7.3 Setup


In order to set up the Vari-Flex jaws for testing, it is first necessary to adjust the height of the tablet shelf on jaw 2. This is done by turning the tablet shelf locking screw counter-clockwise until the shelf moves freely in its guide track. Now place the capsule shaped product horizontally on its side and adjust the shelf up or down so the horizontal center of the product is parallel to the threaded shaft. Then lock the tablet shelf into position by tightening the tablet shelf holding screw with a clockwise rotation.

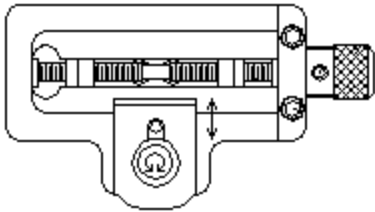

Next, it is necessary to adjust the spacing of the tablet pressure points on Jaw 2. This is done by rotating the thumbscrew located on the side of Jaw 2. This adjustment will move the pressure points closer together or farther apart as needed to test the sample. When making this adjustment, it is recommended that the pressure points be spaced as far apart as possible, while still maintaining contact with the flat surface of the sample.

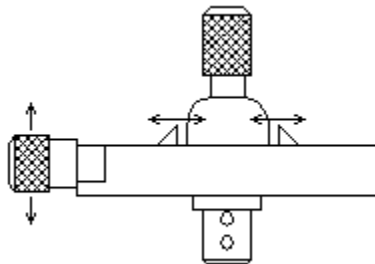


Before putting the machine into operation, verify that the sample is oriented properly, as shown in the illustration to the right. This shows the pressure points behind the sample as dotted lines. The point where the force will be applied by the moving jaw, Jaw 1, is marked by the "X"

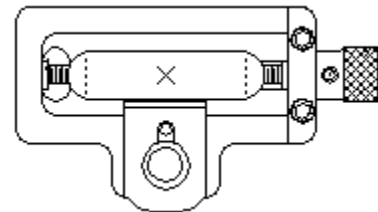

7.4 Operation

After the Vari-Flex jaws have been setup correctly, it is possible to begin testing. The standard operating procedure for the hardness tester may be used at this point. Be sure to watch that the break occurs in the center of the sample as shown in the illustration below. If the sample does not break as indicated or a part of the sample chips off, you may need to change the setup of Jaw 2, so that it makes better contact with the sample. If needed, any faulty sample break can be cleared from the statistical report by pressing the CLEAR key (4) on the keypad of the hardness tester. For more information on the function of the CLEAR key (4) please refer to the Operating Controls section of the manual.

Jaw 1 - Moving jaw (Motor side)

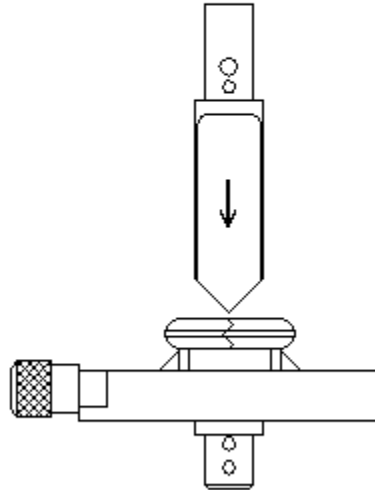

Jaw 2 - Stationary jaw (Load Cell side)

8. TECHNICAL INFORMATION

8.1 Units Conversion List

The following list shows the relationship between the various units that the hardness tester is capable of measuring in:

1 Kilo Pond	=	9.807 Newtons
1 Kilo Pond	=	1.4 Strong Cobbs
1 Strong Cobb	=	0.714 Kilo Ponds
1 Strong Cobb	=	7.005 Newtons
1 Newton	=	0.102 Kilo Ponds
1 Newton	=	0.143 Strong Cobbs
1 Newton	=	2.2048 Lbs.

8.2 RS-232 Serial Port

The connector for the RS-232 serial port of the hardness tester is a 25 pin male Sub-D connector.

8.2.1 Interface

Baud Rate	2400
Parity	Even
Data Bits	7
Stop Bits	1
Handshake	Hardware RTS/CTS

8.2.2 Pin Assignments

Pin	Signal	Description
1	Protective Ground	
2	TXD	Transmit Data
3	RXD	Receive Data
4	RTS	Request to Send
5	CTS	Clear to Send
6	DSR	Data Set Ready
7	GND	Signal Ground
20	DTR	Data Terminal Ready

8.2.3 Data Output

All information printed on the built in printer is automatically sent to the serial port.

8.2.4 Data Input

The following is a list of commands that can be used to control the hardness tester via the serial port:

Command	Function
S	Start Test - same as pressing START key
A	Stop Test - same as pressing STOP key
P	Printer On - turn the printer on
N	Printer Off - turn the printer off

9. OPTIONAL AND SPARE PARTS LIST

Part No.	Description
04-168-0001	CALIBRATION KIT - 98.066 N
04-168-0002	CALIBRATION UPGRADE KIT TO BRING THE 80 N KIT UP TO 98.066 N
15-156-0002A	MEMBRANE COVER FOR KEYBOARD
15-156-0002N	SWITCH FOR KEYBOARD
15-156-0010	HT-500 LOAD CELL
15-156-0014	ELECTROSTATIC PRINTER
15-156-0016	HT-300 LOAD CELL
15-156-0019	LIMIT SWITCH FOR MOTOR HOUSING
15-156-0020	ON/OFF ROCKER SWITCH
15-156-0023	STANDARD JAW SET
15-156-0047	ISOLATION COVER
15-156-0048	PEDESTAL
15-156-0049	PAPER CUTTER FOR PRINTER

Part No.	Description
15-156-0050	PRINthead FOR PRINTER
15-156-0051	POWER CABLE
15-156-0052	PAPER SPINDLE FOR PRINTER
15-156-0053	ELECTROSTATIC PRINTER PAPER
15-156-0055	VARI-FLEX JAWS
15-156-0056	AMPULE HOLDER
15-1561032	KEYBOARD (NEW STYLE)
15-1561033	MAIN BOARD